


CRUTCH HOLDER S3

Click on images for larger format



CRUTCH HOLDER S3

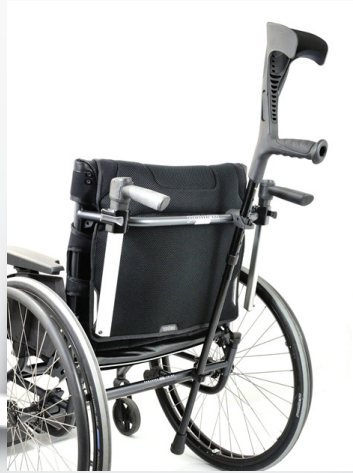

CRUTCH HOLDER S3



LINK TO
DEMO VIDEO

CRUTCH HOLDER S3 - A NEW ACCESSORY [link](#)

Panthera is now launching a crutch holder under its own brand.

The holder consists of two parts that can be easily attached to the crutch. The upper part hooks onto the backrest tubing and the lower attaches to the chassis tubing with a magnet.

The advantage of this crutch holder is that there is no extra weight on the wheelchair and the crutch is very easy to attach and detach from the wheelchair. The crutch is also very firmly attached to the wheelchair.

It is suitable for crutches with a

diameter of 19 to 22 mm (upper attachment) and 16 to 19 mm (lower attachment)

Material: Reinforced polyamide with a recessed magnet

Weight: 46g / crutch

art.no. 9900301

Crutch holder S3