

Which attachment for which model

A basic rule is that wheelchairs with a rigid rear axle, i.e. the rear axle is not adjustable forward or backwards, the attachment is mounted directly onto the rear axle.

NOTE! Except for U3 Light width 45 cm, where (MX2-43P) is used.

On all other models an assembly- kit must be used together with parts of the original attachment.


Model	Attachment (description)	Art.No.	Comment
X	Attachment Smartdrive Panthera U3 Light and X	MX2-44S – 24" MX2-45S – 25" MX2-46S – 26"	Directly onto the rear axle
U3 Light / Light L	Attachment Smartdrive Panthera U3 light och X	MX2-44S – 24" MX2-45S – 25" MX2-46S – 26"	Directly onto the rear axle NOTE! For U3 Light 45 cm you need assembly-kit MX2-43P
S3 / S3 Short	Attachment Smartdrive Panthera U-3/S-3	(MX2-43P)	Complete assembly-kit, see Fig E below
S3 Short- Low	Attachment Smartdrive Panthera U-3/S-3	(MX2-43P)	Complete assembly-kit, see Fig E below
S3 Large	Attachment Smartdrive Panthera U-3/S-3	(MX2-43P)	Complete assembly-kit, see Fig E below
S3 Short Abd	Attachment Smartdrive Panthera U-3/S-3	(MX2-43P)	Complete assembly-kit, see Fig E below
S3 Swing /Swing kort	Attachment Smartdrive Panthera U-3/S-3	(MX2-43P)	Complete assembly-kit, see Fig E below
Extended range			
U3 Light Y Front	Attachment Smartdrive Panthera U3 light and X	MX2-44S – 24" MX2-45S – 25" MX2-46S – 26"	Directly onto the rear axle NOTE! For U3 Light 45 cm you need assembly-kit MX2-43P
U3 Light Y Front 90°	Attachment Smartdrive Panthera U3 light och X	MX2-44S – 24" MX2-45S – 25" MX2-46S – 26"	Directly onto the rear axle NOTE! For U3 Light 45 cm you need assembly-kit MX2-43P
U3 Light Y Front High	Attachment Smartdrive Panthera U3 light och X	MX2-44S – 24" MX2-45S – 25" MX2-46S – 26"	Directly onto the rear axle NOTE! For U3 Light 45 cm you need assembly-kit MX2-43P
U3 Y Front	Attachment Smartdrive Panthera U-3/S-3	(MX2-43P)	Complete assembly-kit, see Fig E below
U3 Y Front Short	Attachment Smartdrive Panthera U-3/S-3	(MX2-43P)	Complete assembly-kit, see Fig E below
U3 Front Low	Attachment Smartdrive Panthera U-3/S-3	(MX2-43P)	Complete assembly-kit, see Fig E below
S3 Swing Low 0°	Attachment Smartdrive Panthera U-3/S-3	(MX2-43P)	Complete assembly-kit, see Fig E below
S3Swing Long	Attachment Smartdrive Panthera U-3/S-3	(MX2-43P)	Complete assembly-kit, see Fig E below

S3 0°	Attachment Smartdrive Panthera U-3/S-3	(MX2-43P)	Complete assembly-kit, see Fig E below
S3 Long	Attachment Smartdrive Panthera U-3/S-3	(MX2-43P)	Complete assembly-kit, see Fig E below

Here is a link to an assembly instruction:

[montering smartdrive faste panthera.pdf](#)

A complete kit contains:

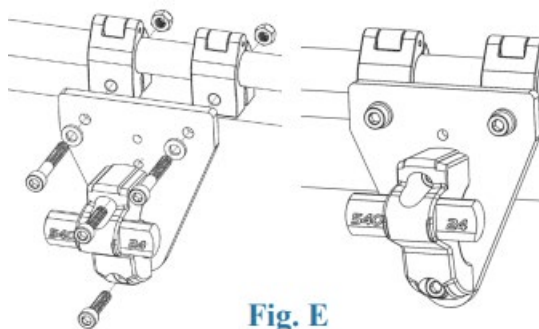

2 p hinge clamp

1p mounting plate

3 p "half-moon" attachment (24", 25", 26")

4p bolts

4p nuts

An attachment for mounting directly onto the rear axle contains:

22" (501) wheel / tyre	Use 24" clamp	
24" (540) wheel / tyre	Use 24" clamp	
25" (559) wheel / tyre	Use 25" clamp	
26" (590) wheel / tyre	Use 26" clamp	

Here is a link to a list over other wheelchair brands/ models compatible/ approved with Smart Drive:

[Wheelchair Compatibilty \(permobil.com\)](#)